

cardiolifecconnect

BECAUSE IT MATTERS  
**SMARTPHONE APP**  
**GPS TRACKING**  
**24/7 AED SERVICE**



# CARDIAC SAFETY

Nihon Kohden's cardiolifeconnect solution provides full-service AED management, including complete monitoring of the cardiolife AED-3100 and other connected AEDs. Not only do you gain insight into the AED's current status, but also its exact location (GPS) and whether it has been used to rescue a patient.

The entire platform is customizable to individual requirements in terms of notifications and integration.

## 24/7 SURVEILLANCE AND READINESS

Cardiolifeconnect is used for 24/7 surveillance of the AED. Nihon Kohden provides monitoring of the AED to verify that the device is operational, such as the shelf life of the disposable pads or the battery level.

You can customize whether you want to receive alerts for daily and monthly self-tests, or when specific actions occur, such as when the AED is turned on or when the AED leaves a certain location.

Cardiolifeconnect allows you to select different notifications.

You can also set whether you want to receive email or SMS notifications about these self-selected events.




E-mail



SMS




Messenger

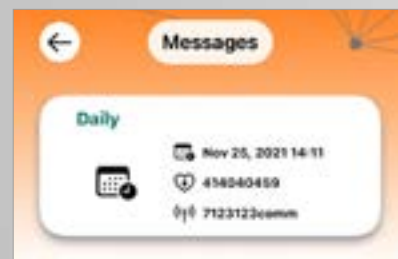

# DASHBOARD AND APPS

Cardiolifeconnect's technology consists of a communication module and a web-accessible dashboard that provides a bird's-eye view of all AEDs, allowing the operator the ability to monitor the service level and readiness of all deployed units at all times. This single, unified view can also be accessed via the CardiLink® apps for iOS and Android smartphones.

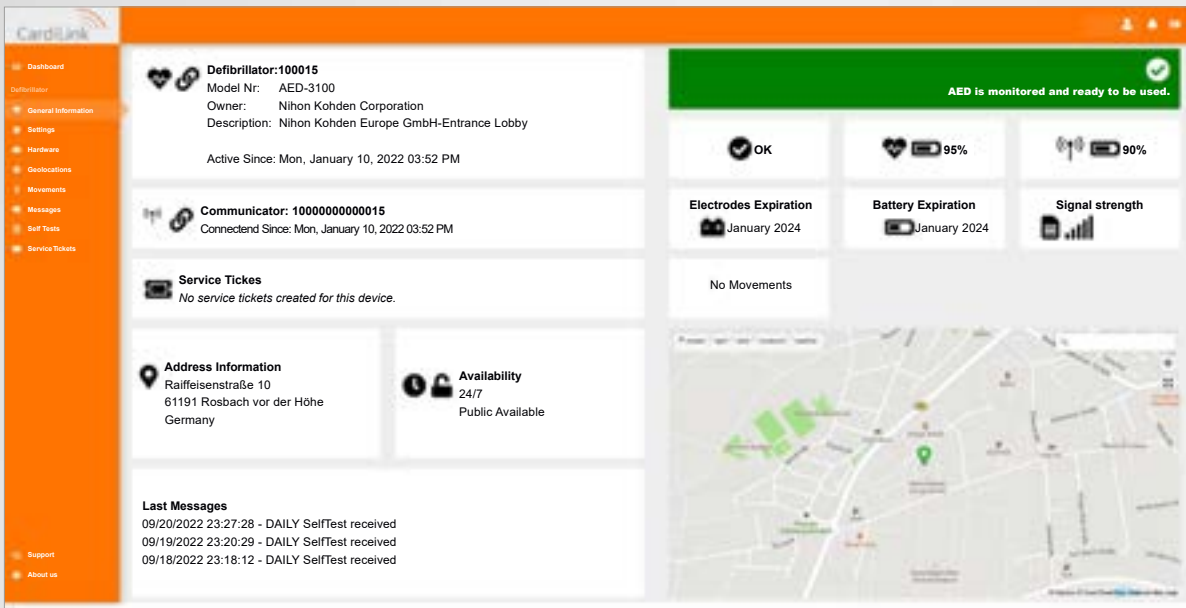

The app intuitively facilitates the data visualisation. Registration is made the easiest way possible via the CardiLink® 2.0 app, available for both Android and iOS devices.





# GPS TRACKING AND GEO FENCING

Cardiolifeconnect enables real GPS tracking. Exact information about the AED's location is displayed on a digital map.



Motion activity details are available in the movement list. The GEO fencing function triggers an alarm if the AED leaves a certain defined radius.

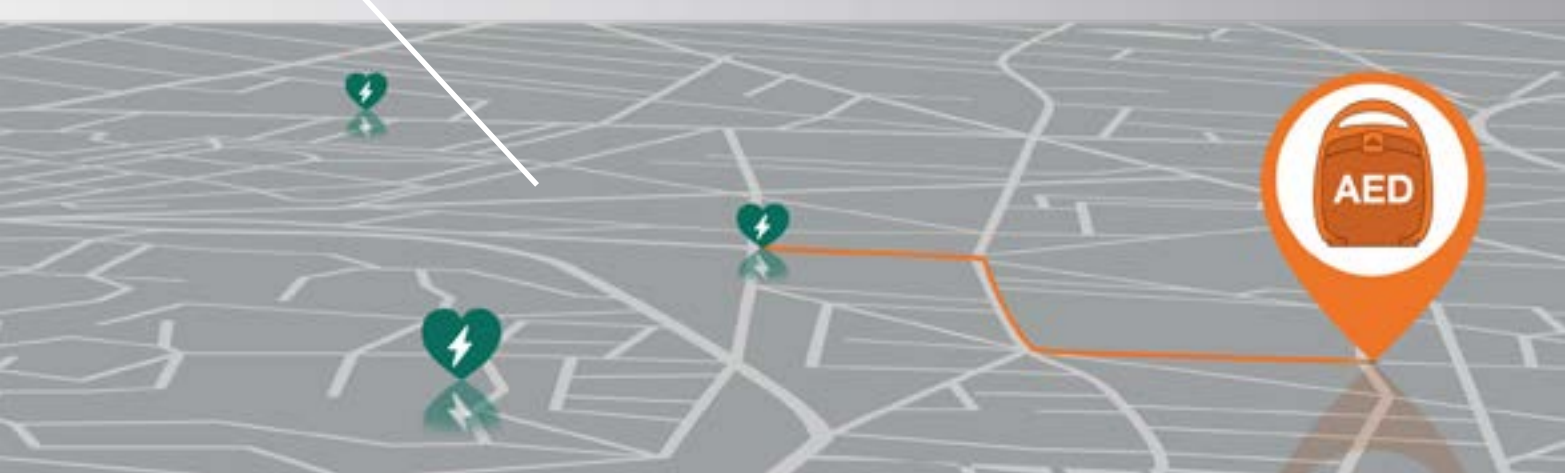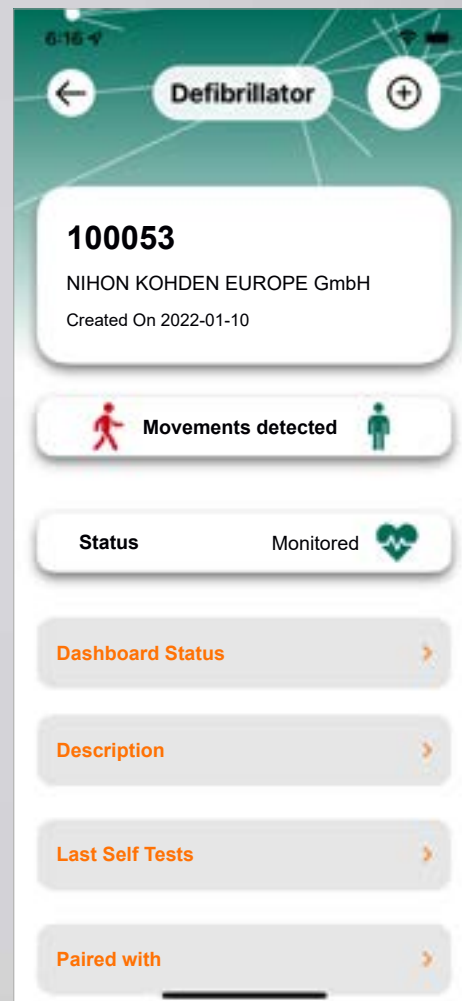




# INTERFACE AND FLEET MANAGEMENT

Due to the interface cardiolifeconnect offers a unique integration which is able to transfer data to mobile first responder applications.

Our fleet management capabilities connect the sometimes disparate technologies deployed by health and safety device manufacturers on one networked platform – mitigating complications arising from local device deployments.

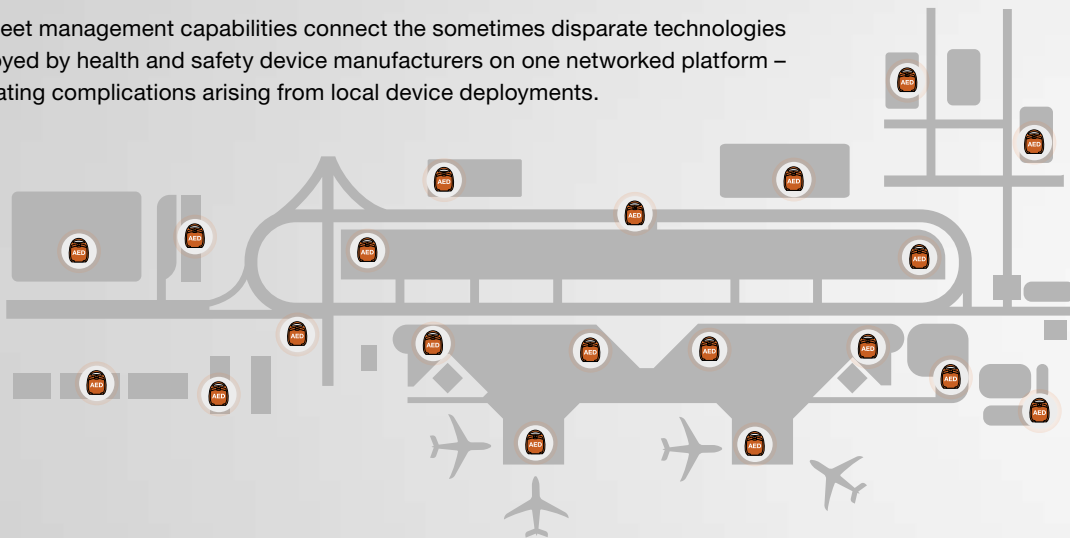

## PRODUCT COMPATIBILITY

This solution offers full hardware compatibility. Nihon Kohden's current cardiolife AED-3100 defibrillator can be equipped with the cardiolifeconnect solution.

cardiolife **AED**  
AED-3100





## Improving Healthcare with Advanced Technology

Since its foundation back in 1951, Nihon Kohden's mission has been to improve the quality of life with advanced technology. We provide solutions for diagnosis, critical care, clinical information, and in vitro diagnostics – and we are dedicated to collaborate with you to meet the challenges of healthcare today and tomorrow.

Visit [www.nihonkohden.com](http://www.nihonkohden.com) to find out more.

#### **NIHON KOHDEN UNITED KINGDOM Ltd.**

Unit 3, Heyworth Business Park  
Peasmarsh, Guildford, Surrey, GU3 1AF, United Kingdom  
Phone: +44 1483 331328,  
Internet: <https://eu.nihonkohden.com>, E-mail: [mail@nihonkohden.co.uk](mailto:mail@nihonkohden.co.uk)

#### **NIHON KOHDEN EUROPE GmbH**

Raiffeisenstrasse 10, 61191 Rosbach, Germany  
Phone: +49 6003 827 0, Fax: +49 6003 827 599  
Internet: <https://eu.nihonkohden.com>, E-mail: [info@nke.de](mailto:info@nke.de)

CardiLink GmbH is the technology partner for  
Nihon Kohden's cardiollifeconnect solution.

#### **CardiLink GmbH**

Gebhardtstraße 28, 6. OG OST, 90762 Fürth, Germany  
Phone: +49 911 733 050 44  
Internet: [www.cardi-link.com](http://www.cardi-link.com)



This document may be revised or replaced by  
Nihon Kohden at any time without notice.