

# Power-PRO TL

powered ambulance cot



# Power-PRO TL

## powered ambulance cot

Our Power-PRO TL powered ambulance cot utilises a battery-powered hydraulic system effectively raising and lowering a cot at the touch of a button.


Carers experience frequent spinal loading due to repetitive motions such as lifting, lowering, carrying and bending. The Power-PRO TL was designed to reduce spinal loading, resulting in reduced injuries, reduced lost and modified workdays and increased recruitment and retention.\*

Each cot lift poses a risk for caregiver injury, but we have the **power** to reduce that risk.

### Key features

- Pneumatically Assisted Backrest. Significantly reduces operator effort and increases patient comfort.
- Cot Fastener System Welded stainless steel construction for durability. Available in floor plate and floor track mounting configurations.
- Wheel-Lock. Power-PRO comes standard with dual wheel locks which help keep cot in place and allows for safer patient transfers.
- 318kg (700lb) capacity.

## Options

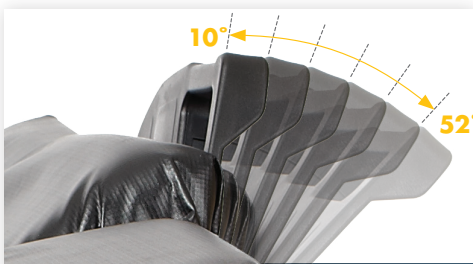

**XPS** expandable patient surface

#### Seven locking positions


Between 10 and 52 degrees. Engineered for durability and cleanability, with a 6061-T6 aluminum core and over-mold design.



**SMRT** power system

#### Automatic charging

SMRT charger automatically evaluates the condition of each SMRT pak and cannot overcharge or overheat.

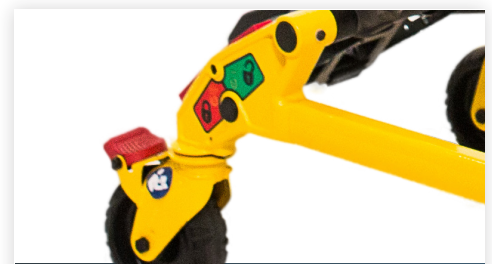

**Steer-Lock** system

Wheel can swivel or lock straight-line position, maximising manoeuvrability.

### Retractable head section

Retractable head section provides 360-degree mobility in any height position.

### XPS

Adjustable with seven locking positions and includes a wider mattress that reduces transfer gap and designed with patient comfort in mind.

### Shock, Flat Leg, or Optional Knee Gatch Positioning

Knee gatch provides patient comfort and greater lift clearance.

### Corner Poles (optional)

Maneuverability in different situations.

### U.S. Patented X-Frame

Will not hot drop.

### Steer-Lock

Maximises cot mobility control by locking the head end casters into a straight position preventing cot drift and enhancing its turning precision.

### SMRT Battery

Zero memory charging solution.

### Hydraulic Lift System

Battery-powered hydraulic lift system with manual back-up for operator confidence.

## XPR restraints

Customise your cot with exclusive, crash-rated patient restraints. Easy to use and quick to clean, we upgraded our legacy design with a polyurethane coating to increase your confidence on scene. XPR restraints are equipped with a polyurethane-coated webbing designed to support cleanability.



### Safety

Enhance your patients' and medics' safety with a proven, crash-rated solution.<sup>1</sup>

Protect the underlying webbing from damage and moisture absorption with an effective polyurethane coating.

### Cleanability

Boost your cleanliness standards to create a more hygienic care environment.

Reduce effort on routine disinfecting protocols.<sup>2</sup>

### Ease of use

Eliminate the learning curve using the same attachment points on Stryker's powered and manual ambulance cots.

Increase productivity and save time with easy-to-read placement indicators.

## Specifications

<b>Model number</b>	6550	
<b>Overall length</b>	206 cm	81 in
<b>Length</b>		
Standard	194.3 cm	76.5 in
Minimum	139.7 cm	55 in
<b>Overall width</b>	58 cm	23 in
<b>Maximum weight capacity</b>	318 kg	700 lb
<b>Height range (to litter top)</b>		
High	106 cm <sup>3</sup>	41.5 in
Low	36 cm	14 in
Weight	66 kg <sup>4</sup>	145 lb
Backrest articulation	0° - 75°	

### ProCare Services

Every day, you count on your medical equipment to perform at its best. With ProCare Services, our highly trained technicians help ensure your equipment is ready to perform when needed. Our goal is to make life easier for your staff and clinicians with service plans designed to help you get the most out of your investment. Choose the service plan that meets your needs.

### ProCare Protect

A complete service solution designed to maintain your equipment with annual preventive maintenance while protecting your investment with an unlimited number of repairs, including parts, labour and associated travel.

### ProCare Prevent

Supports your preventive maintenance initiatives with an annual inspection of your equipment, including parts labour and associated travel.

### Flex Financial Program

Our financial programs provide a range of smart alternatives designed to fit your organisation's needs. We offer flexibility beyond a cash purchase with payment structures customisable to meet budgetary needs and help to build long-term financial stability. Contact your Stryker representative for more information.

Program availability will vary depending on country and product portfolio, please check with your local sales representative for full details and availability.

1. Meets SAE J3027 and BS EN 1789:2020 dynamic crash standards in combination with Power-LOAD and Power-PRO XT, Power-PRO TL or Performance-PRO XT. 2. When power washed within cot system. 3. Height measured from bottom of mattress, at seat section, to ground level. 4. Cot is weighed with one battery pack, without mattress and restraints.

Stryker reserves the right to change specifications without notice.

\*Evaluation of Medical Cot Design Considering the Biomechanical Impact on Emergency Response Personnel – T.K. Fredericks, S.E. Butt, K.S. Harms, J.D. Burns XXVth Annual Occupational Ergonomics and Safety Conference 2013.

This document is intended solely for the use of healthcare professionals. A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area. Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Flex Financial, Performance-LOAD, Power-LOAD, Performance-PRO, Power-PRO, ProCare, XPR and Stryker. All other trademarks are trademarks of their respective owners or holders. The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

