



ProVu™ Lite

Video Laryngoscopes

Get a fresh perspective

Thoughtfully designed for consistent visualisation.



Reusable controller with locking single use low-profile blades.



Designed for both routine practice and difficult airways without compromise.



High quality image display featuring auto-focus and anti-fog.

ProVu™ Lite Video Laryngoscope

Designed for both routine practice and difficult airways without compromise.

ProVu Lite extends familiar laryngoscopy practice with consistent visualisation.

The landscape orientation of the display provides clear anatomical and airway device visualisation, enabling minute adjustments to positioning, ensuring an accurate and successful placement.

Integrated visualisation

ProVu Lite combines familiar laryngoscopy technique with a purpose-built video architecture that expands where and how visualisation can be used.

Ready in an instant

ProVu Lite has an intuitive “plug and go” interface. Simply connect the display and controller to wake. The responsive, text-free display allows you to focus on the camera output.

70° Adjustable display tilt

When space is restricted, the ProVu Display can be tilted, optimising your visual placement of the laryngoscope blade without changing your technique.

Anti-Fog design

The internally mounted Anti-Fog design reduces the need for warm up time and allows users to intubate quickly with a clear view.

High performance LED

A fully encased, bright LED light source targets the glottic area to deliver a clear visual of the airway. Suitable for both video intubation and direct laryngoscopy techniques.

Specifications

	Controller	Blade
Sterility	Non-Sterile	Sterile
Number Of Uses / Reprocessing Cycles	2,500	Single Use
Connection	Micro USB	USB-C
Ingress Protection	IP44	-
LED Illuminance	-	>500 Lux



Start every case with a new lens

By integrating a new camera into every disposable blade, ProVu Lite provides consistent image quality from case to case, eliminating long-term optical degradation and variability seen with reusable camera systems.

This approach supports clinician confidence during everyday intubations while maintaining readiness for more challenging airways.

Reusable 3.5" or 8" display

Reusable controller

Single use handle and blade with camera & LEDs



Get everyone involved

Achieve the optimal screen viewing angle. Connect ProVu Lite to a ProVu 8" Display and mount on the ProVu Display Pole Stand and Universal Articulating Arm, improving visibility for clinicians and their team, especially during training and demonstrations.

Storage for ProVu devices in the basket and charging capabilities on the pole.



The ProVu Lite Controller is a reusable, lightweight component, universally connecting the single use blade to the display, enabling seamless video transmission and intuitive device control.



+ 70° Adjustable screen tilt

When space is restricted, the connected ProVu Display or amplifier cable can be tilted, optimising your visual placement of the laryngoscope blade without changing your technique.

+ Easy cleaning and intermediate-level disinfection

Smooth surfaces enable rapid reprocessing using common hospital wipes and chemicals.

+ Handle locking mechanism

Eliminates accidental ejection during use. Requires an intentional, two-finger, press to release or eject the blade following the intubation procedure.

+ Engineered for consistent, reliable performance

Rigorously tested for 2,500 reprocessing cycles, enduring repeated insertions, hinge rotations and display connections to ensure durability.



Bespoke blades & cameras

Each ProVu Lite video laryngoscope blade features precisely angled, optimally positioned camera placement to deliver a clear, reliable view.

Rather than relying on a single reusable camera inserted into disposable sheaths, ProVu Lite features a **dedicated camera** integrated into every blade. This ensures that the optics and illumination are **tailored** to the specific blade geometry.

With every blade, clinicians receive a new, size-specific camera aligned to provide the highest-quality image, enhancing airway visualisation and maximising first-pass success.

Our low-profile **Macintosh** blades are ideal for routine intubation, featuring the familiar shape as our successful and reliable BriteBlade Pro blades, enhanced with the advantages of video technology to aid intubation.

Each **Hyperangulated** blade provide a superior view of the glottic opening without having to align the oral, pharyngeal and laryngeal axis, reducing the need for lifting pressure and potential for trauma to the patient. Ideal for difficult, anterior airways and patients with restricted neck movement.

Each blade incorporates:

- + Optimised camera geometry
- + Sealed optics
- + Integrated LED illumination
- + Reinforced structure
- + New camera in each single use blade
- + Blind-spot reduction technology
- + Low profile design
- + Blind-spot reduction technology
- + Clear glottic exposure
- + 60° Field of view
- + Anti-fog technology
- + Anti-glare



ProVu™ Display Pole Stand + Universal Articulating Arm

Your all-in-one
intubation workstation

Experience effortless airway visualisation with a complete, ready-to-use intubation workstation designed to elevate every procedure.

Easily position your display to personalised, precise angles for rapid decision-making and seamless team communication, with multiple adjustment fixtures.

Keep your ProVu devices organised and within reach using the integrated storage basket, and stay powered at all times with convenient built-in charging.



Key Features

- Adapt to your environment
- Achieve optimal viewing angle
- Mount twin arms to multiply video outputs capabilities
- Locking wheels



ProVu Reusable 3.5" Display

Each reusable 3.5" ProVu Display battery lasts for 90 minutes of continuous use.



ProVu Reusable 8" Display

Battery lasts for up to 3.5 hours of use from a fully-charged battery, or can be plugged in for continuous use.


ProVu Articulating Arm
with Adjustable Clamp



ProVu Display
Table Stand




ProVu Universal Reusable
2m Amplifier Extension Cable

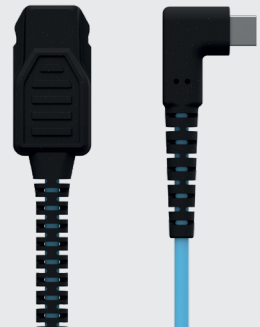

ProVu 3.5" Display
with Charging Dock



Protective Casing
for 3.5" Display

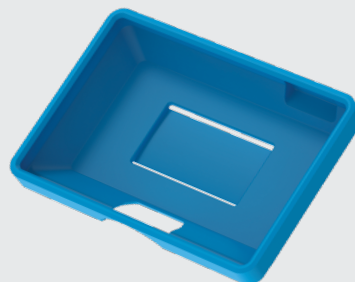

ProVu™ Lite

CODE	DESCRIPTION	QTY
040-11-0000	ProVu Lite Video Laryngoscope – Reusable Controller	1
040-12-0130	ProVu Lite Video Laryngoscope – Macintosh 3 Handle and Blade	10
040-12-0140	ProVu Lite Video Laryngoscope – Macintosh 4 Handle and Blade	10
040-12-2130	ProVu Lite Video Laryngoscope – Hyperangulated 3 Handle and Blade	10
040-12-2135	ProVu Lite Video Laryngoscope – Hyperangulated 3.5 Handle and Blade	10
040-12-2140	ProVu Lite Video Laryngoscope – Hyperangulated 4 Handle and Blade	10

ProVu™ Displays & Accessories

CODE	DESCRIPTION	QTY
040-07-0035	ProVu 3.5" Display + Charging Dock	1
040-07-0080	ProVu 8.0" Display + Charging Cable + 2m Amplifier Cable	1
040-07-0020	ProVu Universal Reusable 2m Amplifier Extension Cable	1
038-990-101	ProVu Display Pole Stand + Universal Articulating Arm	1
038-990-108	ProVu Articulating Arm with Adjustable Clamp	1
038-990-102	ProVu Display Table Stand	1
038-990-017	ProVu Reusable Silicone Protective Casing for Universal 3.5" Display	2

ProVu™ Rigid Stylets

CODE	DESCRIPTION	QTY
040-08-2024	ProVu Rigid Stylet – Medium	10
040-08-2030	ProVu Rigid Stylet – Large	10



Flexicare (Group) Limited

Cynon Valley Business Park, Mountain Ash, CF45 4ER, UK

+44 (0)1443 474647 • enquiries@flexicare.com

www.flexicare.com

