

iPAD **NFK200** **Defibrillator**

This is a new model of the iPAD series public access AED. The **NFK200** is a simple to use, intuitive defibrillator allowing for quick and convenient pads access.



**AFFORDABLE
EFFECTIVE**

KEY FEATURES

- Easy access to pads from front panel
- Device and consumables status LED indicator
- CPR metronome, voice guidance, and graphic

INSTRUCTIONS

- Data transfer by USB

SAFETY

- Automatic internal discharge
- Daily / weekly / monthly self test
- Adult / paediatric mode change button

TECHNOLOGY

- Semi automated e cube biphasic defibrillation
- Combined adult / paediatric pads
- CPR step detection indicator for more effective CPR



www.welmedical.com

NFK200 Defibrillator

IPAD

The Important Key Features

NFK200

Device and consumables status LED indicator

Adult/ Paediatric switch mode change button



Combined adult / paediatric pads

CPR metronome & voice guidance

Easy access to pads from front panel

Graphic Instructions

Data transfer by USB

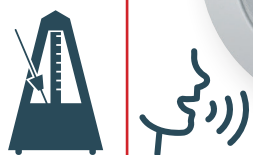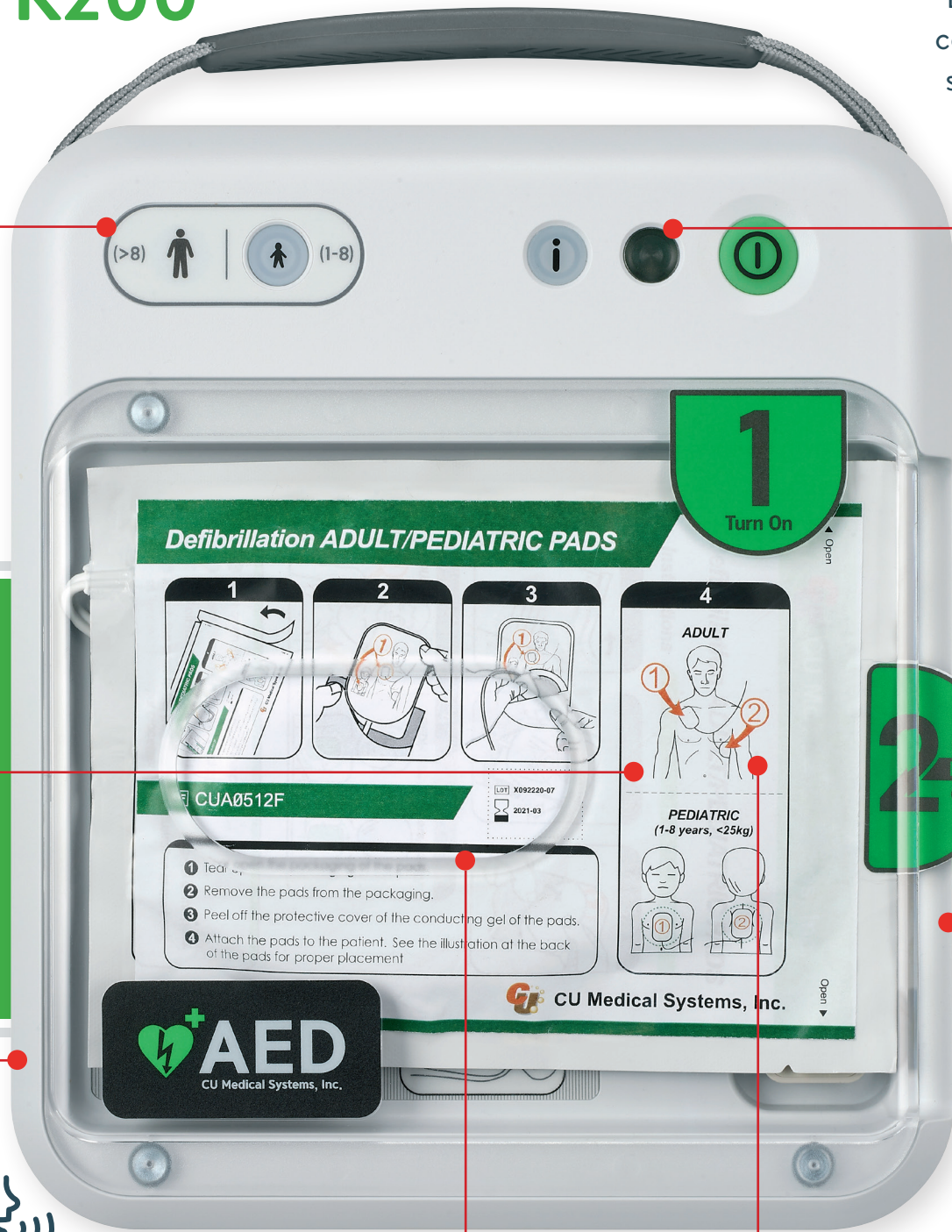

An affordable compact system

Features Include:

- Rapid adult / paediatric change button
- Rugged casing for use in all conditions, and handle for easy carry
- Supports CPR and resuscitation of various patients in diverse locations

The **NFK200** is a simple, compact, highly effective addition to the iPAD series. It can easily and effectively be used by anyone.

SUITABLE
FOR
THE
WORKPLACE



EASY
TO USE

The NFK 200 Defibrillator is available in the **UK only** from **WEL Medical** and **selected business partners**.

This will be a popular option for those seeking an affordable and effective defibrillator.

Easy access to pads from front panel



Contact us now if you need more information or if you would like to place an order:

T 01252 344007

Physical

Dimensions	260mm x 220mm x 70mm (Width x Length x Height)
Weight	2.2kg (Including the battery pack and pads)

Environmental

Operating Conditions	Temperature: 0°C ~ 50°C (32°F ~ 122°F) Humidity: 5% to 95% (non-condensing)
Storage Conditions	Temperature: 0°C ~ 50°C (32°F ~ 122°F) Humidity: 5% to 95% (non-condensing)
Transport Conditions	Temperature: -20°C ~ 60°C (-4°F ~ 140°F) Humidity: 5% to 95% (non-condensing)
Altitude	0 to 4,572 m (0 to 15,000 ft.)
Drop	Withstands 0.7 meter drop to any edge, corner or surface
Vibration	Operating: Meets MIL-STD-810G Fig.514.6E-1, Standby: Meets MIL-STD-810G Fig.514.6E-2
Sealing	IEC 60529:2013 IP42
ESD	Meets IEC 61000-4-2:2008
EMI (Radiated)	Meets IEC 60601-1-2 limits, method EN 55011:2016+A1:2017, Group 1, Class B
EMI (Immunity)	Meets IEC 60601-1-2 limits, method EN 61000-4-3:2006+A2:2010 Level 3 (10V/m 80MHz to 2.5GHz)

Defibrillator

Operating Mode	Semi-automated
Waveform	E-cube biphasic (Truncated exponential type)
Output Energy	150J at 50 Ω load for adults 50J at 50 Ω load for children
Charge Control	Controlled by an automated patient analysis system
Charging Time	Within 3 seconds from when the voice instruction, "An electric shock is needed" is issued.
Time from initiation of rhythm analysis	10 seconds with a new battery (even after the delivery of 15 discharges at 150J) 12 seconds with a new battery (even after the delivery of 15 discharges at 200J)
Time from Power on to readiness for discharge	25 seconds with a new-fully charged battery (even after the delivery of 15 discharges at maximum energy)
Time from CPR to Shock	At least 6 seconds from the completion of CPR to shock delivery
Disarm	Patient's heart rhythm changes to non-shockable rhythm The SHOCK button is not pressed within 15 seconds

ECG Acquisition

Acquired ECG Lead	Lead II
Frequency Response	1 Hz to 30 Hz

ECG Analysis System

Impedance Range	25Ω to 175Ω
Shockable Rhythms	Ventricular Fibrillation or Fast Ventricular Tachycardia
Sensitivity and Specificity	Satisfies AAMI DF80

WEL MEDICAL

1 CHANCERYGATE WAY · FARNBOROUGH · HAMPSHIRE · UK · GU14 8FF

Controls, Indicators, and Prompts

Controls	Power Button Shock Button Adult/Pediatric Selection Button
Indicators	Do-Not-Touch-Patient Pads Patch Position Indicators Status LED Indicator CPR Step and Detection Indicator Shock Button
Speaker	Provides voice prompts
Sound Level	75dB, 1m above speaker
Beeper	Provides various audible indications
Low Battery Indicator	Shown on the Status LED, announced via voice instruction

Self-Test

Automatic	Power On Self-Test, Run-time Self-Test Daily, Weekly, and Monthly Self-Test
Manual	Battery Pack Insertion Test

Battery Pack Long Life

Battery Type	12V DC, 4.2Ah LiMnO ₂ , Disposable
Capacity	150 Joules - At least 200 Shocks for a new battery or 8 hours of operating time
Standby Life	At least 5 years from the date of manufacture
Temperature Ranges	Operating Temperature: 0°C ~ 50°C (32°F ~ 122°F) Storage Temperature: -20°C ~ 60°C (-4°F ~ 140°F)

Adult / Paediatric Defibrillation Pads

Surface Area	85cm ²
Cable Length	120cm
Shelf Life	At least 36 months from the date of manufacture

Data Storage and Transfer

USB	External memory. Data may be copied from the internal memory to the USB.
File System	FAT32
Internal Memory Data Capacity	5 individual treatments, up to 3 hours per treatment



T 01252 344007

orders@welmedical.com

DEFIBRILLATORS | DEFIBRILLATOR DEVICES | DEFIBRILLATOR TRAINING | CONSUMABLES

www.welmedical.com