



Manual Power Assist Bikes

Empulse F55Product Overview





8.5" Wheel

14" Wheel





Empulse F55



Power add-ons are designed to enhance the mobility of user.

Current options are:

- Permanently fitted onto the chair and don't mirror the energy patterns of users
- Require core strength to fit onto the chair
- Designed solely for performance and suit only very active clients.

Empulse F55 has been designed to work for semi-active to active clients.

Blending the benefits of add-on wheels and push devices which enhance clients mobility with the technology and function of a powered bike.

[Removable, easy mounting and intuitive to use]





Key Features



Simple Mounting

- Only mounted on the chair when required mirroring energy patterns
- Intuitive to mount for clients with limited core strength
- Fits on both folding and rigid chairs with fixed or swing-away front frames



Everyday Performance

- Compact and lightweight for indoor usage and transportation
- Designed for everyday use optimum speed, range and performance
- Two wheel sizes 8.5" solid or 14" pneumatic wheel



Suitable for Everybody

- Intuitive controls and handling for the best driving experience
- Lots of built-in adjustment to suit a wide range of clients
- A boxed product for simple ordering and set up





Simple Mounting

Easy Mounting

- Simple to mount to any chair allowing clients to fit the Empulse F55 when required throughout the day. Meaning clients don't have to propel extra weight, saving energy throughout the day
- Fits to the front of the chair making it simple for clients to handle versus products
 that mount to the rear of the chair. An optional kick stand is available allowing the
 F55 to stand on its own for clients with poor trunk control making it even easier to
 dock.









Easy Mounting Continued

- Tapered frame receivers mean it's more intuitive to fit the subframe arms into the receivers than other powered add-on bikes
- The F55 is light to **lift and lock** into position making it simpler for clients with limited strength to work with; this feature also lifts the castors from ground meaning they don't get caught in cracks in the pavement.






Simple Mounting

Universal Mounting - Frame Types

Suitable for both rigid and folding chairs with either a fixed or swing-away front frame.

Designed to suit semi-active to active users, the F55 ensures the handlebar is in the correct position irrelevant of the chair frame type. This is achieved by offering a swingaway stem which has more depth than the fixed front option.



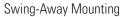



Fixed Front Mounting

Universal Mounting - Fits All Chairs

Available with 5 clamp options to fit to a wide range of chairs:

- Round clamp for 19.5, 23, 25, 28.6 & 30mm diameter tubes (Neon² or Argon²)
- Oval clamps for QUICKIE chairs (Xenon² FF, Nitrum or Helium)
- Xenon² Swing Away clamps
- Oval clamps for Progeo (Joker, Ego & Tekna)
- Ki Mobility Rogue clamp (35mm diameter)





Empulse

Everyday Performance - Indoors & Transport

1. Indoor Performance

Designed for everyday use in everyday tasks, the F55 allows the bike to sit closer to the chair to reduce the overall length. Meaning its easy to navigate tight homes, shops or malls.

An optional 8.5" wheel provides a smaller package for the best indoor performance.



2. Transport

Ideal for transport due to its lightweight and compact design.

From just 11kg it can be easily lifted and stored in the car boot.

The subframe can be locked up or removed for an even smaller and lightweight package.

Certified for aircraft transport the F55 can also be taken abroad on holidays.









Everyday Performance - Outdoors







Suitable for Everybody - Controls

Intuitive Controls



Why *complicate* things?

We designed the F55 with ease of use in mind, catering for the first time user rather than our more advanced clients

1. LCD Display

Standard with a clear LCD display. Which is simply to function and clearly shows the speed display and other mode options.

2. Forward & Reverse

Changing direction is easy, simply push the 'direction' switch up or down to easily move from forwards to reverse.

LED Light

Comes as standard with a bright built in LED light simply hold down the + button for 2 seconds to activate making those early morning and late night journeys a lot brighter.

3. E-Brake

As well as the traditional disk brake we have introduced an electric brake for a softer, smoother more controlled stop to the journey.

4. Cruise Control

Built in to the base price is a cruise control feature allowing clients to relax on those longer journeys by maintaining a safe speed.







Suitable for Everybody - Adjustable

Designed to fit a wide range of clients and chairs the Empulse F55 is a boxed product which is easy and quick for assessments and handovers; everything needed is in the box or adjustable on the bike.

Empulse F55 Adjustments

- 1. Frame width adjustment 26 48cm Includes 4 centre bar lengths
- 2. Handle bar stem height 15cm quick release
- **3.** Handle bar angle 70° quick release
- **4. Subframe Length** 30cm
- 5. Clamp height and 360° clamp angle adjustment








Easy Ordering

Included in the Empulse F55 base price:

- 8.5" or 14" drive wheel
- 6, 10 or 15kph motors
- LCD Display
- LED Light
- Forward & reverse
- Cruise control
- Electronic brake
- 4 frame width bars
- 4 colour accents red, blue, orange and black
- 8.1Ah Lithium Battery
- Charger
- Chair clamps

Simple to order with only 5 options to select choose from:

- Drive wheel size 8.5" or 14"
- 2. Front frame type fixed or swing-away
- Clamp type QUICKIE or third party
- 4. Speed Setting 6kph, 10kph or 15kph
- 5. Accessories Kick stand, extra battery or charger

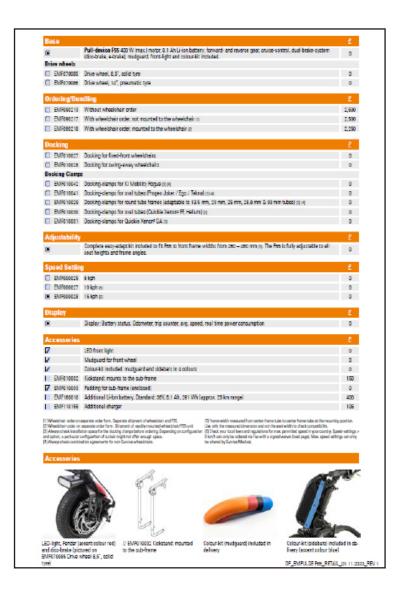





Commercials

1. Empulse F55 only

When the client is fitting the F55 to their current day chair:

- RRP £2,500
- Dealer NET £1,500
- Dealer Margin £1,000



2. Empulse F55 bundling with QUICKIE chair

Available when bundling with a NEW QUICKIE chair allowing the user to save £250 on the Empulse F55

- RRP £2,250
- Dealer NET £1,350
- Dealer Margin £900





WEIGHTS & USE

Maximum user weight

110 kg

F55 weight

From 13 kg

Transport weight

From 9 kg

Removable parts

Sub frame and battery

Crash tested

Yes

Transport without user

Remove or lock off sub frame so unit is flat and store securely in vehicle

Intended use environment

For people unable to walk short or long distances, with limited mobility both inside and out



Sunrise Medical Limited

Thorns Road, Brierley Hill, West Midlands, DY5 2LD

Tel.: +44 (0) 845 605 66 88 Fax: +44 (0) 845 605 66 89

E-Mail: enquiries@sunmed.co.uk www.sunrisemedical.co.uk

In line with our policy of continual product improvement, Sunrise Medical reserves the right to make changes at any time, and without notice.